


NOTES:

- ① JUNCTION BOXES: J1--NORMAL POWER, J2--EMER. POWER, J3--TELEPHONE.
- ② GAS PIPING: TYPE L COPPERTUBING. OXYGEN & AIR--1/2" I.D., VACUUM--3/4" I.D.
- ③ COLOR TO MATCH 'SURF' OR 'PLATINUM GREY'.
- ④ CEILING HEIGHT: 7'6" MIN TO 9'6" MAX. OTHER HEIGHTS ARE SUBJECT TO AN ADDTL CHARGE. CONTACT YOUR SALES REPRESENTATIVE FOR DETAILS.
- ⑤ MAXIMUM CONDUIT SIZE FOR DOUBLE-SIDED UNITS IS 3/4".

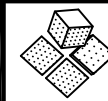

MODULAR 500 EAST BRITTON ROAD
SERVICES COMPANY OKLAHOMA CITY, OK. 73114

PROJECT



REVISED

PROFILE® UNIT
 SINGLE SIDED, SEMI-RECESSED
 DUAL GAS POD
 MODEL 7011-A

SCALE: = 1'


DATE _____ DRAWN BY _____ DRAWING #
 PB-109

SERVICES TO BE REVERSED
 ON OPPOSITE HAND UNITS

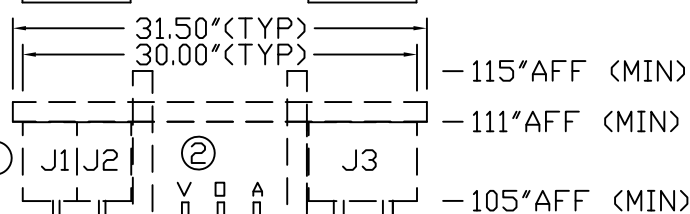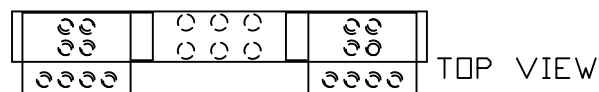

MEDIUM ACUITY


-- REQ'D: -- AS SHOWN
 -- OPPOSITE

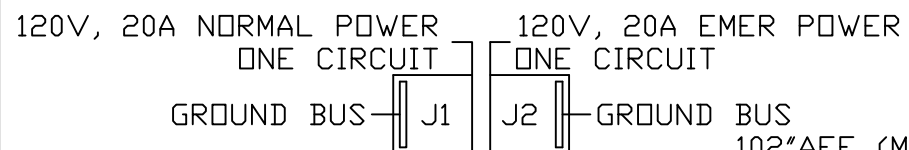


102" AFF (MAX)

ACCESSIBLE
 CEILING LINE

90" AFF (MIN)

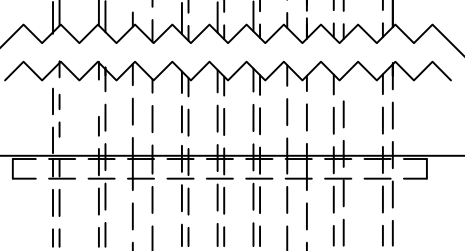

VTF

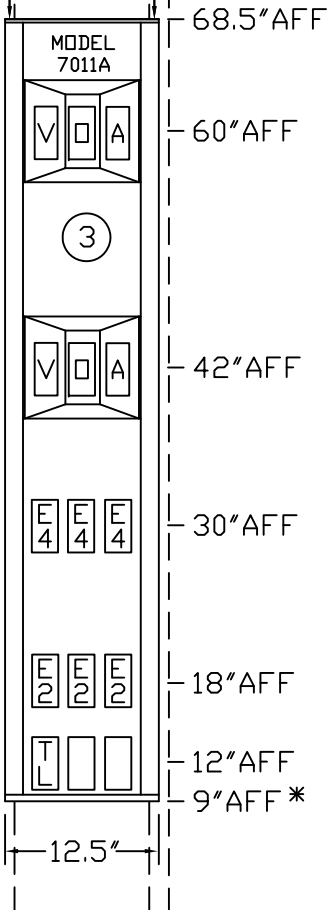

STRUCTURAL
 WALL

3.625"



2.5"

1"

3.25"

SIDE VIEW

CIRCUIT 1
 EMERGENCY POWER

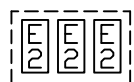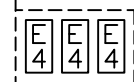

CIRCUIT 1
 NORMAL POWER

WIRING DIAGRAM
 GROUND WIRING TO BE #10
 GREEN STRANDED COPPER

FLOOR LINE

*AFF (ABOVE FINISHED FLOOR)

