


- NOTES:
1. FORM UNIT IS SEMI-RECESSED. SERVICES ARE TOP FED. UNIT IS BUILT IN TWO VERTICAL SECTIONS.
 2. THE ACTUAL LOCATIONS OF THE JUNCTION BOXES, CONDUIT TERMINATION POINTS AND MEDICAL GAS TERMINATION POINTS ARE SUBJECT TO CHANGE.
 3. SINGLE SIDED UNITS INSTALL IN A 3-5/8" PARTITION; BACK TO BACK UNITS INSTALL IN A 6" PARTITION. ALL UNITS SET FOR A 5/8" GYPSUM BOARD FINISH.
 4. INSET SURFACE PROVIDED ABOVE AND BELOW PANELS FOR GYPSUM BOARD TO BE PROVIDED AND INSTALLED BY OTHERS.
 5. INTER CONNECTIONS MAY BE REQUIRED BETWEEN MODULES.
 6. SERVICES TO BE REVERSED ON OPPOSITE HAND UNITS

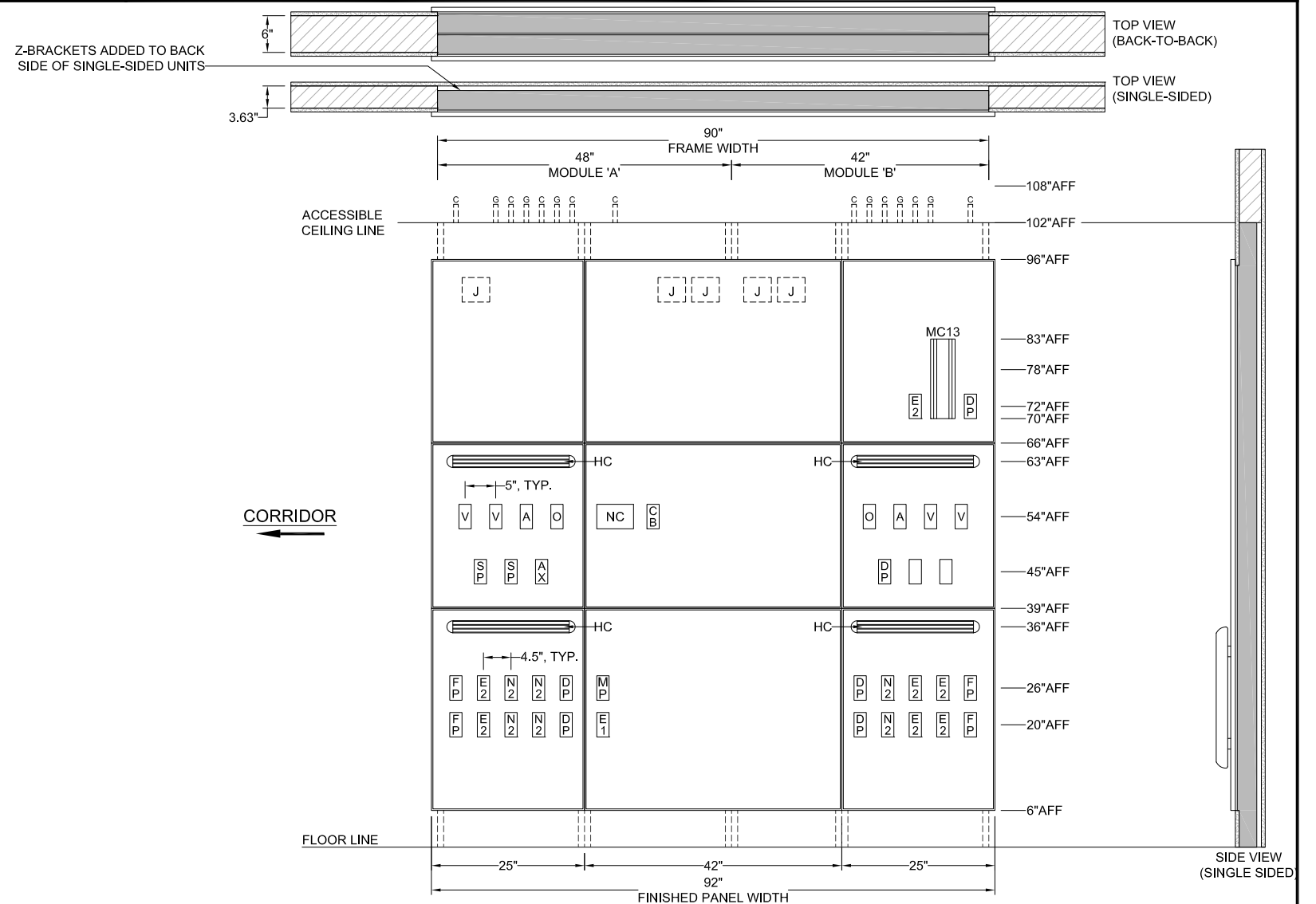

GENERAL SPECIFICATIONS:						LEGEND:				ACCESSORIES NOT SHOWN (PER UNIT)			
ELECTRICAL		PLUMBING		FINISHES		REQ	SYM	DESCRIPTION	REQ	SYM	DESCRIPTION	REQ	DESCRIPTION
NORMAL POWER JUNCTION BOX:	J	COPPER TUBING TYPE:	L	LAMINATE: L1		2	O	OXYGEN OUTLET	6	DP	DATA OUTLET PROVISION		
EMERGENCY POWER JUNCTION BOX:	J	OXYGEN SIZE (I.D.):	1/2"	TYPE:	LAMINATE BRAND:	2	A	AIR OUTLET	4	FP	FUTURE PROVISION		
LIGHTING JUNCTION BOX:	J	MEDICAL AIR SIZE (I.D.):	1/2"	-	-	4	V	VACUUM OUTLET	1	MP	MULTI-PIN BED RECEPT. PROVISION		
LV CONDUIT STUBS:	C	VACUUM SIZE (I.D.):	3/4"	BACK LIT LED RESIN PANEL: R1 OR LAMINATE: L2		6	N2	DUPLEX RECEPTACLE, 8300 IVORY 20A, 125V, H.G., NORMAL POWER	1	MC13	MONITOR CHANNEL ONLY, 13" LENGTH		
CONDUIT SIZE; LINE VOLTAGE:	1/2"			TYPE:	LAMINATE BRAND:	1	E1	SINGLE RECEPTACLE, 8310 RED, 20A, 125V, H.G., EMER. POWER	4	HC	HORIZONTAL ACCESSORY CHANNEL		
CONDUIT SIZE; LOW VOLTAGE:	3/4"			-	-	7	E2	DUPLEX RECEPTACLE, 8300 RED 20A, 125V, H.G., EMER. POWER					
CONDUIT CONNECTORS:	STD. SET SCREW					1	NC	NURSE CALL PROVISION					
WIRE TYPE:	THHN, STRANDED			ELECTRICAL COVER PLATES: STAINLESS OR NYLON TO MATCH DEVICE		1	CB	CODE BLUE PROVISION					
WIRE SIZE:	12 AWG			HORIZONTAL CHANNELS: ANODIZED ALUMINUM		1	AX	AUXILIARY ALARM PROVISION					
GROUND WIRE COLOR:	GREEN			UNIT TRIM: ANODIZED ALUMINUM		2	SP	SWITCH PROVISION					

<p>500 EAST BRITTON ROAD OKLAHOMA CITY, OK 73114 www.modularservices.com</p>	TOTAL QTY:	X	FORM® 7500 SERIES SEMI-RECESSED	
	AS SHOWN:	X	DRAWN BY:	XXX
	OPPOSITE:	X	DRAWING #:	4XXXX
PROJECT:	HOSPITAL PROJECT CITY, STATE	DATE:	XX/XX/XX	SALES ORDER #: XXXX XXXX